



**MEGA
MEDICAL**

www.megamedical.com.au

**OXYGEN CONCENTRATORS +
MEDICAL OXYGEN CONSUMABLES**

Medical Oxygen Consumables	1
Medical Oxygen Masks	1
Medical Oxygen Tubing	3
Accessories and Parts	3
Oxygen Concentrators	4
Portable Concentrators	4
Stationary Concentrators	8

Introduction

Mega Medical is Australia's leading supplier of medical gas equipment, with more than 20 years of experience in the industry.

Australian-owned and operated, Mega Medical proudly serves all Australian states and territories health services and supplies medical equipment nationwide. We provide our customers with access to high-quality products, backed by the excellence of our customer sales and support division.



**MEGA
MEDICAL**

www.megamedical.com.au

Medical Oxygen Consumables

MEDICAL OXYGEN MASKS

ADULT NASAL CANNULA WITH 2M TUBING

Brand: Salter Labs

SKU: LR3233

Description

The nasal cannula is used to deliver supplemental oxygen to a patient in need of respiratory help. Also referred to as 'Nasal Prongs' these are used to administer oxygen to the patient via the nose to increase efficiency in gas delivery and provide a more comfortable and convenient delivery method.

Patient Notes: This product is designed for single patient use and is not intended to be reprocessed. Reusing, soaking, rinsing or sterilising this device with chemicals may leave harmful residues and may impede the functions of this device. Reusing single patient use device may also increase the risk of cross contamination. It is strongly advised that the product be used according to trained medical personnel's guidance and the instructions for use before the marked expiration date.

Specifications

- Length: 2m
- Colour: Clear
- Material: Medical PVC

Features

- Soft, transparent, odour free vinyl to ensure the best fit and user comfort.
- Flexible oxygen tubing prevents occlusion when tubing is kinked.

Standard Contents

- Soft tip adult nasal cannula
- Patient cautionary notes



Related Products

50

Adult Nasal Cannula with 2m Tubing - Box 50
SKU: LR3233-Box50



Child Nasal Cannula with 2m tubing
SKU: LR3232

50

Child Nasal Cannula with 2m tubing - Box 50
SKU: LR3232-Box50



OXYGEN MASK WITH 2M TUBING – ADULT

Brand: Salter Labs

SKU: LR3645

Description


Adult oxygen therapy mask or 'Hudson Mask' as it is commonly referred to, is the ideal solution in patient oxygen therapy. This oxygen therapy mask is manufactured by Salter Labs and is one of the premium masks on the market. Moulded from superior odour-free vinyl, this mask is made for ultimate function and comfort. The mask features a thick elastic strap and an adjustable aluminium nose clip to ensure the most secure and comfortable fit for the patient. The 2 metre translucent tubing will connect to all medical oxygen therapy units with an oxygen outlet nipple.

Features

- Soft, transparent, odour-free vinyl
- Aluminium fixer and elastic strap to ensure the best fit and user comfort

Related Products

50

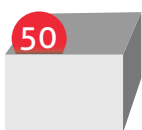

Oxygen Mask with 2m Tubing - Adult - Box 50


SKU: LR3645 - Box50



Oxygen Mask with 2m Tubing – Child

SKU: LR3644

50

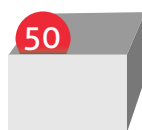

Oxygen Mask with 2m Tubing - Child - Box 50

SKU: LR3644-Box50



Child High-Oxygen Concentration Mask

SKU: LR3674



Adult Oxygen Concentration Non-Rebreathing Mask

SKU: LR3675

MEDICAL OXYGEN TUBING



Oxygen Tubing 10m
SKU: LR31610



Oxygen Tubing – 2m
SKU: LR2002



Oxygen Tubing Connector
SKU: BOS1215

ACCESSORIES AND PARTS



Oxygen Humidifier Bottle
Disposable
SKU: BOS7600



Oxygen Tubing Y Connector /
Flow Splitter
SKU: PAAPM88409

Oxygen Concentrators

PORTABLE CONCENTRATORS



ECLIPSE 5 PORTABLE OXYGEN CONCENTRATOR

Brand: Caire

SKU: 6900AU-SEQ

Description

The Eclipse 5 combines portability with clinical efficiency to treat a wide variety of your needs, no matter the range of your activity levels – giving you the freedom to travel and enjoy life on the go. Suitable for use 24/7 use, this all-in-one oxygen therapy device can provide continuous flow from 0.5 to 3 LPM (Litres Per Minute) and pulse doses up to a setting of 9.

Continuous Flow Settings (LPM):

0.5, 1.0, 1.5, 2.0, 2.5, 3.0

Pulse Dose Settings (mL):

16, 24, 32, 40, 48, 56, 64, 72, 80, 88, 96, 128, 160, 192

Features

PATIENTS

- Run on AC, DC or battery power
- In the car—all flow rates available and battery will recharge ≤ 2 LPM
- 2 hours battery at 2 LPM with single battery
- Over 5 hours of battery duration at 2 pulse
- Removable rolling cart with adjustable handle
- FAA approved for use on commercial airline

Specifications

DIMENSIONS

- Height: 49cm
- Width: 31cm
- Depth: 18cm
- Weight: 8.3 kg with battery

OTHER

- 2 hours of battery at 2 LPM with single battery
- Removable rolling cart with adjustable handle
- Runs on wall outlet (AC), motor vehicle (DC), or rechargeable battery power

Includes

- Transportable oxygen concentrator
- Power cartridge (battery)
- AC power supply
- DC power supply
- Removable cart
- Accessory bag
- Nasal cannula
- Spare intake filter

Related Products



Eclipse 5 Protective Cover
SKU: 5052-SEQ



Eclipse 5 - Power Cartridge (Battery)
SKU: 7082-SEQ



CAIRE FREESTYLE COMFORT CONCENTRATOR – 16 CELL BATTERY

Brand: Caire

SKU: AS200-104

Description

Expand your comfort zone with this innovative and powerful device. The Freestyle Comfort is lightweight, has a curved design that hugs your body, so you will hardly know it's there. It has the power to let you stay out all day long with greater flexibility to increase flow rates as your prescription changes. Patented smart oxygen delivery features ensure you get the oxygen you need, when you need it, keeping you comfortable and worry-free by adjusting to your automatic breathing patterns and level of activity.

Please note: The therapeutic use of Medical Oxygen requires a prescription, trained clinical assessment and proper instruction on use. Medical Equipment & Gases Australia Pty Ltd requires a copy of the patients prescription before the product can be dispatched

Specifications

- Flow settings: 1-5 pulse dose only
- Dimensions with 16-cell battery pack: 28.8 x 18.5 x 7.9 cm
- Weight (with 16-cell battery pack): 2.6 kg
- Sound level: 39.9 dBa @ Setting 2
- Battery duration with 16-cell battery pack at setting 2 = 8 hours
- Warranty: 3 year warranty on FreeStyle Comfort (includes sieve beds); 1 year warranty on battery and power supplies

Features

- Ergonomic design for optimal comfort
- Easy-to-Read Display – Large fonts, multi-color LCD display and glow-in-the-dark keypad.
- Meets FAA requirements for commercial airline travel
- Up to 16 hours battery life with the 16-cell battery pack
- Delivers 1-5 pulse flow settings
- Maximum output of up to 1050 mL per Minute
- 2.3 kgs with a small, compact frame

Standard Contents

- Portable Concentrator
- Carry Bag
- 16 Cell Battery Pack
- AC/DC Universal Power Supply
- Shoulder Strap
- UltraSense Technology
- Warranty: 3 years unit including sieve beds. 1 year battery & power supplies

Pulse Settings

- 1 = 210 mL/min
- 2 = 420 mL/min
- 3 = 630 mL/min
- 4 = 840 mL/min
- 5 = 1050 mL/min

Related Products



Freestyle Comfort - 8 Cell Replacement Battery

SKU: BT034-1S



Caire Freestyle Comfort Universal Power Supply

SKU: ASPW023-3

Variations



Caire Freestyle Comfort Concentrator – 8 cell battery
SKU: AS200-4



ZEN-O – PORTABLE OXYGEN CONCENTRATOR

Brand: GCE

SKU: CIGRS-00502

Description

Zen-O™ is a best in class portable oxygen concentrator with pulse and continuous oxygen delivery modes for active patients requiring oxygen therapy. Zen-O™ features breath sensitivity and Rate Response Therapy™ (RRT) that automatically adjusts to the patient's demands.

Zen-O™ delivers up to 2 litres of oxygen per minute in response to the patient's need. Unlike other devices that deliver a fixed amount of oxygen, Zen-O™ automatically increases the amount of oxygen delivered if a patient's breath rate rises.

Includes: Zen-O concentrator, 2 x 12 cell lithium-ion battery, AC power supply, car charger, power cord, pull cart, carry bag, user manual, quick start guide

Specifications

- Size (W×D×H): 212 mm × 168 mm × 313 mm (8.3" × 6.6" × 12.3")
- Weight: 4.66 kg (10.27 lbs) with one 12 cell battery
- Battery duration: Up to 4 hours with a single battery or 8 hours with 2 batteries at 18 breaths per minute
- Power requirements: AC adaptor: 100-240V AC(+/- 10%) 50-60 Hz in, 24V DC, 6.25 out; DC adaptor: 11.5 – 16V DC in, 19V, 7.9A out
- Purity: 87% – 96% at all settings
- Maximum oxygen discharge pressure: 20.5 psi
- Inspiratory trigger sensitivity: -0.12cm/H2O
- Humidity range: 5% to 93% ± 2% non-condensing
- Operating altitude: 0' to 13,000' (0m to 4,000m)
- Operating temperature: 5°C (41°F) and 40°C (104°F)
- Storage temperature: -20°C (-4°F) and 60°C (140°F)
- Setting: Adjustable in 0.5 increments from 1.0 to 6.0 in pulse mode and from 0.5 to 2.0 in continuous mode (LPM)
- Noise level: 38 dB(A)
- Alarm types: Low oxygen purity, no breath detected, low battery and service required
- Warranty – Device: 3 years or 15,000 hrs
- Warranty – Sieve beds: 2 years or 15,000 hrs
- Warranty – Battery and accessories: 1 year

Features:

- **DUAL FLOW MODE** Zen-O™ offers patients the best of both worlds. Patients can alternate between continuous flow and pulse mode oxygen therapy.
- **SIMPLE AND EASY TO USE** Zen-O™ is designed with patients in mind, it is simple to use with intuitive button operation and LCD display.

Related Products



Zen-O battery 12 cell
SKU: RS-00501



SIMPLYGO PORTABLE OXYGEN CONCENTRATOR

Brand: Philips Respironics

SKU: RES1069058

Description

SimplyGo is Philips first oxygen concentrator that offers the benefits of a continuous flow mode in a lightweight, easy-to-carry design of 4.5kg. No need to compromise between an effective oxygen therapy output and the usual weight of the device.

SimplyGo offers a solution that fulfills the needs of most oxygen therapy patients. With SimplyGo, patients have the flexibility to be away for the day, for work, for a weekend or long holiday.

Standard

Pulse mode bolus size:

- 1 = 12ml
- 1.5 = 18ml,
- 2 = 24ml
- 2.5 = 30ml
- 3 = 36ml
- 3.5 = 42ml
- 4 = 48ml
- 4.5 = 5.4ml
- 5 = 60ml
- 5.5 = 66ml
- 6 = 72ml

Features

- Control panel display
- Power – Turns SimplyGo on and off
- Alarm Silence – Turns the audible alarm off
- Mode – Selects one of the three possible operating modes
- Display screen – Shows information about the operating status of the device
- Plus (+) / minus (-) – Increased / decreases the displayed setting



Respironics SimplyGo Mini
- Standard Battery
SKU: RES1113604



Respironics SimplyGo Mini
- Extended Battery
SKU: RES1113605

Related Products



SimplyGo Mini Standard
Lithium Ion Battery
SKU: RES1116816



SimplyGo Mini Extended
Lithium Ion Battery
SKU: RES1116817



VISIONAIRE 5L OXYGEN CONCENTRATOR

Brand: Caire

SKU: AIR-Vis5

The VisionAire 5 oxygen concentrator is one of the market leading 5L units available. Coming from the reputable oxygen machine manufacturer 'Caire', their brand of units are market leading and the VisionAire 5 is no different. Built with greater portability in mind the 5-litre unit, which is the smaller version of the NewLife units, can deliver great reliability at a fraction of the cost and weight. Lightweight, quiet and power efficient, this stationary oxygen concentrator is the must-have for your home oxygen delivery needs.

The unit is simple to use with a clear and easy to use operating switch to turn the unit on and off, as well as a clear and reliable flowmeter and control knob to safely and easily adjust the oxygen flow rate. The unit comes with a power cord that meets Australian standards and is as simple as plug in and switch on!

Recommended accessories are the 'LR3233 – Adult Nasal Cannula with 2m Tubing'



Specifications

- Flow Rates: 1.0 – 5.0 LPM
- Oxygen Concentration: 90 % (+5.5 / -3 %)
- Dimensions: (52.8 x 35.8 x 29.2 cm)
- Weight: 30 lbs (13.6 kg)
- Sound Level: 45 dB(A)
- Operating Environment: 41 to 105°F (4 to 40°C), 15-95% humidity
- Storage Environment: -4 to 140°F (-20 to 60°C), 15-95% humidity
- Power Consumption: 290 W (All Flow Rates)
- Alarms: Power Failure, High Pressure, Low Pressure, No Flow, Low O2 Concentration (Optional)
- Warranty: 3 years
- Maintenance Schedule: None
- Max Tubing Length: 50 ft (15.2 m)

Features

- Lightweight 13.6kg frame that fits easily into any home.
- Power-efficient, drawing less than 300 watts at all continuous flow settings minimising expenses.
- A large humidifier bottle stand that allows for easy attachment of both re-fillable and pre-filled humidifier bottles.
- User-interface contains a back-lined flowmeter for viewing accuracy and proper setting.
- Equipped with a resettable circuit breaker, a dual-function hour meter and well-defined system alarms.

Standard Contents

- 1x VisionAire 5L Oxygen Concentrator
- 3 Years Limited Warranty
- 1x User Manual



NEWLIFE INTENSITY 10 LITRE DUAL FLOW OXYGEN CONCENTRATOR

Brand: Caire

SKU: AS099-200

HIGH FLOW AND HIGH PRESSURE COMBINED

Combining high flow and high output pressure makes the Intensity the premier high flow stationary concentrator for at home and facility use. This concentrator is also perfect for use in Veterinary Clinics as a portable oxygen supply for sedation.

High-Pressure Applications:

The Intensity has an output pressure of 20 psi (137.9 KPa), more than double any other high flow concentrator.


Specifications

- Flow Rates: 2-10 LPM
- O₂ concentration LPM: 92%+3.5%/-3%, (10 LPM: 90%+5.5%/-3%)
- Dimensions: 69.9 x 41.9 x 36.8 cm
- Output Pressure: 20 psi (138 kPa)
- Power Consumption: 590 W
- Weight: 58 lbs (26.3 kg)
- Sound Level: 55 dB(A)
- Operating Environment: 5°C to 40°C
- Max Tubing: 200 ft (61 m)

Features

- Warranty 3 Years
- Oxygen Therapy
- Dual Flow
- Provide a total of 10 LPM
- Allows two patients to get oxygen from one machine simultaneously
- Perfect for hospitals, long term care facilities
- Maintenance Schedule: Felt Filter – 1 Year Replacement; Intake Filter – Clean Weekly

Related Products

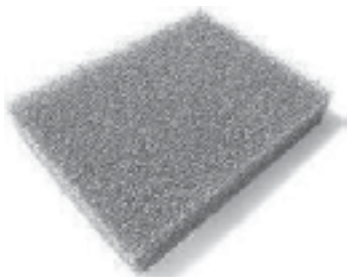

NewLife Intensity Foam Filter - 15PPI
SKU: F1002-1



NewLife Intensity 10 Litre Single Flow Oxygen Concentrator
SKU: AS099-211



EVERFLO STATIONARY OXYGEN CONCENTRATOR

Brand: Philips Respironics

SKU: RES1020010

Description

The EverFlo from Respironics is a unique stationary oxygen concentrator that delivers the features home care providers want and patients deserve. This top of the line unit from Philips is the benchmark in quality home oxygen concentrators. With a robust and compact design, this Everflo can produce a flow rate from 0.5 to 5 Litres per minute all while operating quietly and consistently.

The unit is simple to use with a clear and easy to use operating switch to turn the unit on and off, as well as a clear and reliable flowmeter and control knob to safely and easily adjust the Oxygen flow rate. The unit comes with a power cord that meets Australian standards and is as simple as plug in and switch on!

Recommended accessories are the 'LR3233 – Adult Nasal Cannula with 2m Tubing'

Maintenance

- Air Inlet Filter is to be changed every 24 MONTHS
- The filter door should be cleaned MONTHLY or as environmental conditions warrant.

Standard Contents

- 1x Respironics EverFlo Oxygen Concentrator
- 3 Years Limited Warranty
- 1x EverFlo User Manual

Features

- Lightweight and easy to transport at 14kgs
- Compact, ergonomic design takes up less space
- One of the quietest mainstream oxygen concentrators
- No homecare provider filter change for two years
- Uses less power and produces less heat, so it's more patient-friendly

Specifications

- Input Voltage: 240 V
- Input Frequency: 50 hz
- Average Power Consumption: 350 watts
- Oxygen Concentration (at 5 LPM): 93% (+/-3%)
- Litre Flow: 0.5 to 5 litres per minute
- Weight: 14 kg (31lbs) – (Shipping weight 17kgs)
- Sound Level: 45 dBA typical
- Dimensions: Height 584mm (23')
- Dimensions: Width 381mm (15')
- Dimensions: Diameter 241mm (9.5")
- Required Maintenance for Air Intake Filter: Change every 2 years
- OPI (oxygen percentage indicator) alarm levels: Low Oxygen – 82%, very low oxygen – 70%

Related Products



EverFlo – Air Inlet Filter
SKU: RES1038831



**MEGA
MEDICAL**

www.megamedical.com.au

Medical Equipment & Gases Australia

ABN 76 053 999 205

9 Middleton Road Cromer NSW 2099

1300 881 824 info@megamedical.com.au